

Cub[®] Pediatric Crib

Safe. Access. Ease of use.

When you're using our products that are developed with safety and ease-of-use at the forefront of design, you can feel confident you're giving your patients the best care possible. The Cub Pediatric Crib has a simple, safe design that's warm friendly appearance cultivates a comforting environment for small children and their families.



Swing-open access doors allow convenient patient care without disturbing a sleeping child. The access doors have built in safety features to help ensure a child cannot open it from the inside.



All four siderails can be lowered to provide unobstructed access which allows caregivers to have 360° access.



All activation points are colored red and require two adult size hand movements to initiate. This helps to ensure movements are deliberate and the patient remains in the crib.

Retractable fifth wheel helps improve caregiver mobility.

Retracting backrest helps eliminate head-end gap while raising the head end of the stretcher.

Counterbalance system on siderail for easy raising and lowering.

Optional protective top can be equipped with intermediate support position on each side.


Optional integrated scale system obtains accurate weights with the child in any position.*

Multiple sticker color options allow you to customize the stretcher for identification purposes.



*The scale system allows the reading of a patient's weight in pounds (lb) or kilograms (kg) with a ± 0.5 lb (0.2 kg) at 50 lb or less ± 1.01 lb (0.4 kg) at ≥ 50 lb (22 kg).

Cub Pediatric Crib

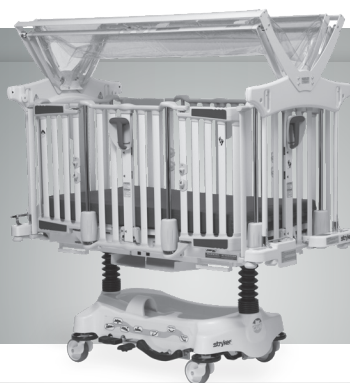

Standard features

- 4" (10 cm) Ultra Comfort support surface
- Manually-operated pneumatic-assist backrest
- Manual foot prop
- Siderail with 9" (23 cm) mandatory stop
- Counterbalance system on siderail for easy raising and lowering
- Oxygen bottle holder
- Access doors at all four corners
- Fully-lowerable side rails provide unobstructed access around entire unit
- Retractable fifth-wheel steering
- 6" (15 cm) Omni Surface™ casters
- Central-locking brake system
- Dual side-mounted foot control hydraulics
- Uni-lower pedal
- Hands-free Trendelenburg and reverse Trendelenburg

Warranty

- One year parts, labor and travel or two years parts

Stryker reserves the right to change specifications without notice.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Cub, Stryker. All other trademarks are trademarks of their respective owners or holder.

3800 E. Centre Avenue
Portage, MI 49002 USA
t: 269 329 2100
toll free: 800 327 0770

patienthandling.stryker.com

Specifications

Model number	FL19H
Overall length and width	
No accessory supports	65.5" x 37" (166 x 94 cm)
With accessory supports	72.25" x 37.5" (184 x 95 cm)
Weight capacity	400 lb. (181 kg)
Safe working load	100 lb. (45 kg)
Height range (to litter top)	
High	40" (102 cm)
Low	33" (84 cm)
Litter positioning	
Manual activated	0° / 25° / 40° / 55°
Pneumatic assist	0 - 50°
Trend./reverse Trend.	±9.5°
Foot section angle	0° / 6° / 12°
Support surface	30.13" x 58.25" (76.5 x 148 cm)
Caster diameter	6" (15 cm)

Optional features

- Integrated scale system (scale system is not CE marked)
- Siderail without 9" (23 cm) mandatory stop
- Accessory supports with roller bumpers (head, foot or both ends)
- Chart shelf
- IV caddy
- Monitor tray (not compatible at foot-end when scale option is selected)
- Oxygen bottle retaining ring
- Siderail pads
- IV poles
 - Permanent folding two- or three-stage
 - Standard removable
- Retractable protective top with or without intermediate supports
- Upright oxygen bottle holder