

stryker

Performance-PRO™ XT

manual ambulance cot



Performance-PRO™ XT

manual ambulance cot

As a medic, the unknown is everywhere. Uncertain situations can arise on every call. We have the power to help prepare for uncertainty. Introducing, our Stryker Performance-PRO XT ambulance cot. Its aluminum construction and 700 lb weight capacity afford a lightweight and versatile manual ambulance cot with enhanced stability.

Performance-PRO XT's three-position adjustable load height accommodates multiple ambulance heights and encourages proper body mechanics when loading and unloading. The optional XPS siderails increase surface area of the cot by a maximum of 38%, offering a solution for patients of a variety of sizes. The retractable head section allows for 360° mobility in any height position for maximum versatility. Duplicated head and foot end lift and grip controls provide ergonomic lifting positions that accommodate caregivers of different sizes.

Versatility. It's a powerful thing

The Performance-PRO XT is our tough, aluminum manual X-frame cot offering. Versatility, proven performance, and increased stability define our Performance-PRO XT cot.




Key features

- Lift-capable safety bar assures handling confidence. Reduces lift height for smaller operators.
- Hand clearance around the litter frame allows for operator safety.
- Adjustable load height industry-exclusive, three position settable load height adjusts to different ambulance heights.
- Built-in pull handle offers maximum operator control during transport and stows out of the way when not in use.



Side Release Handle

Side release handle allows one-person raising and lowering of an unoccupied cot.

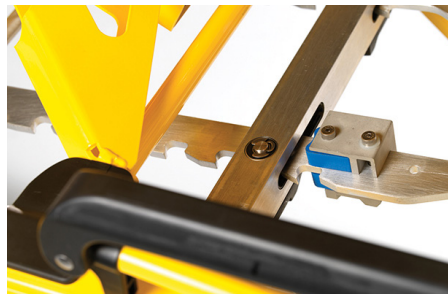

Adjustable Load Height

Three-position-settable load height adjusts to different ambulance heights.



Head and Foot end Lift and Grip Section

Assists in ergonomic lifting.

XPS

Adjustable with seven locking positions and includes a wider mattress that reduces transfer gap and designed with patient comfort in mind.

Pneumatically assisted Backrest

A Stryker innovation. Significantly reduces operator effort and increases patient comfort.

Shock, Flat Leg, or Optional Knee Gatch Positioning


Knee gatch provides patient comfort and greater lift clearance.

X-Restraint Package

Meets SAE J3027 dynamic crash test safety standards.

Retractable Head Section

Retractable head section provides 360-degree mobility in any height position.

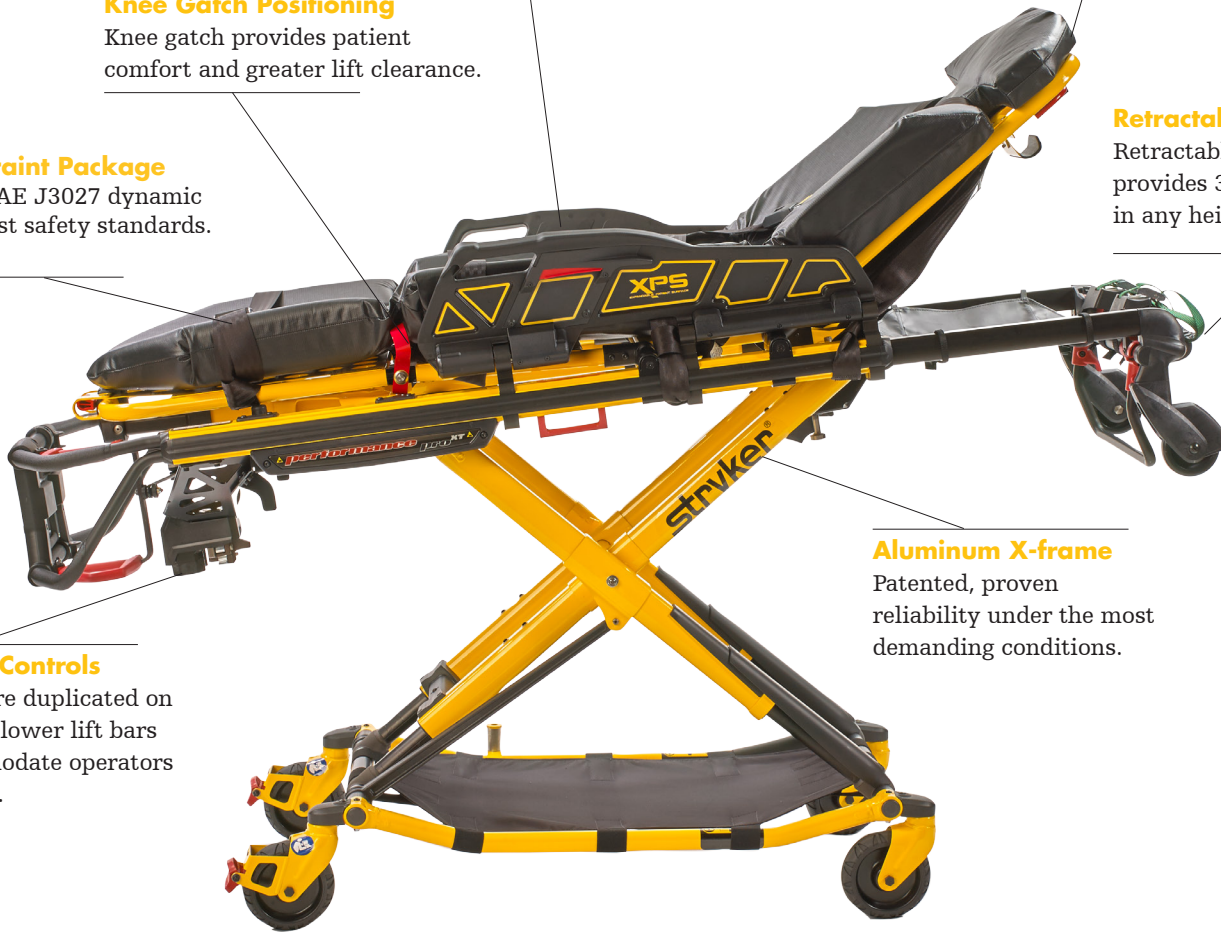

Aluminum X-frame

Patented, proven reliability under the most demanding conditions.

Foot end Controls

Controls are duplicated on upper and lower lift bars to accommodate operators of all sizes.

up to **38%**

increased surface area with XPS

700LB

weight capacity

360°

degrees of mobility in any height

Specifications

Model Number	6086
Height¹	
Position 1	13.8 in (35.1 cm)
Position 2	22 in (55.9 cm)
Position 3	25.8 in (66.5 cm)
Position 4	28.1 in (71.4 cm)
Position 5	31.9 in (81 cm)
Position 6	34.6 in (87.9 cm)
Position 7	37.3 in (94.7 cm)
Position 8	40 in (101.6 cm)
Position 9	42.2 in (107.2 cm)
Load Height Adjustment	
High	Deck heights 34 in (86 cm) and above
Mid	Deck heights between 30 – 34 in (76 – 86 cm)
Low	Deck heights below 30 in (76 cm)
Length	
Standard	80 in (203 cm)
Minimum	64 in (163 cm)
Width	23 in (58 cm)
Weight²	89 lb (40.37 kg)
Wheels	
Diameter	6 in (15 cm)
Width	2 in (5 cm)
Articulation	
Backrest	0 – 73°
Shock Position	+ 15°
Optional Knee Gatch	30°
Maximum Weight Capacity	700 lb (318 kg)
Minimum Operator Required	
Occupied Cot	2
Unoccupied Cot	1
Recommended Fastener System	
Power-LOAD	Model 6390
Floor Mount	Model 6370 or 6377
Wall Mount	Model 6371
Recommended Loading Height³	Up to 34 in (86.4 cm)

Technical Support

Stryker Technical Support comprises a team of professionals available to help with your Performance-LOAD needs. Contact via phone at 1 800 STRYKER or email at TechnicalSupport@Stryker.com.

ProCare Service from Stryker

With maintenance, service and work-flow solutions that proactively maintain your equipment, we help you get every last penny out of your investment. Two-hour priority response callback time and 24- to 72-hour product turnaround help you focus more on your patients and less on equipment upkeep.

Flex Financial Program

Our financial programs provide a range of smart alternatives designed to fit your organization's needs. We offer flexibility beyond a cash purchase with payment structures customizable to meet budgetary needs and help to build long-term financial stability. Specialized finance solutions for municipalities and other government agency businesses are also available.

To learn more about the wide range of financial services, call 1 888 308 3146.

¹Height measured from bottom of mattress, at seat section, to ground level.

²Cot is weighed without mattress and restraints.

³Load height can be set between 27.25 in and 34 in.

Stryker reserves the right to change specifications without notice. In-service video included with every order.

The Performance-PRO XT is designed to conform to the Federal Specification for the Star-of-Life Ambulance KKK-A-1822.

The Performance-PRO XT is designed to be compatible with competitive cot fastener systems.

The yellow and black color scheme is a registered trademark of Stryker Corporation.

Meets dynamic crash standards with Performance-LOAD cot fastening system (SAE J3027 with X-restraints).

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: **Flex Financial, Performance Pro, Power-LOAD, ProCare Stryker**. All other trademarks are trademarks of their respective owners or holder.